


## Operating instructions

### How to treat your pimples and pustules

Pimples and pustules are both treated the same way with CLARO. Each treatment consists of two 6-second applications - 12 flashes per application.

To treat a pimple or pustule, complete these 4 steps:

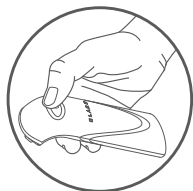

#### STEP 1: Turn On

- Press and Release Button
- Green Light Turns On




#### STEP 2: Get Ready

- Place Firmly on Pimple
- Listen for 2 Beeps
- Keep Firmly on Pimple



#### STEP 3: First Application


- **PRESS BUTTON AGAIN**  
*without removing from skin*
- Wait for 12 Flashes
- Lift the CLARO from your skin

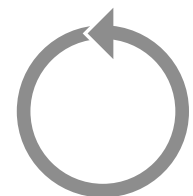

#### STEP 4: Second Application

- Repeat Steps 2 and 3 right away
- You are done with the complete treatment

## Having Trouble? Useful User Hints

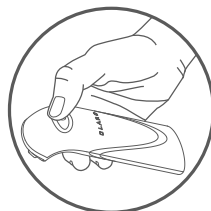

#### HINT 1:

- Press the **Control Button** once quickly and release. The **Light Ring** turns green when the CLARO is ready to use.
- “Quickly” means less than 1 second. If you press the Control Button longer than 1 second, the unit may not turn on.

#### HINT 2:

- Center the **Treatment Window** over the pimple or pustule. Make certain that the **Skin Contact Sensors** are in full contact with your skin. Wait until you hear the 2 beeps.
- Use a mirror to help center the treatment window over the pimple or pustule you are treating. You are ready to begin.
- If skin contact is lost, the treatment will end immediately. If this happens, start over at Step 2.



#### HINT 3:

- To begin the treatment, press and release the **Control Button**. Maintain skin contact at all times. Each application will consist of 12 light flashes and last 6 seconds.
- After an application is completed, the Light Ring will turn red.
- Keep CLARO in contact with skin.
- One complete treatment = Two 6-second applications for a total of 24 flashes.

